



# TAKE TREATMENTS WHEREVER LIFE TAKES YOU!



**Take treatments wherever life takes you**

with the new PARI TREK<sup>®</sup> S Compact Compressor Nebulizer System. With its lightweight design, multi-voltage AC adapter, and handy carrying case, you can take the PARI TREK<sup>®</sup> S with you wherever you go.

Equipped with the PARI LC<sup>®</sup> SPRINT Reusable Nebulizer, the PARI TREK<sup>®</sup> S delivers powerful performance and shorter treatment times.

PARI's smallest and most compact compressor weighs less than a pound, and at approximately 50 dB, the PARI TREK<sup>®</sup> S is one of the quietest compressors on the market.




Ambient transport temperature	-13°F to 158°F (-25°C to 70°C)	Weight	(compressor) 0.8 lbs. (360 gms) (with optional battery) 1.2 lbs. (520 gms)
Humidity	10% to 95%	Dimensions	(compressor) 4.5 x 3.5 x 1.75 inches
Operating time with battery	Approx. 50 minutes	Electrical requirements	(compressor) 100VAC - 240VAC 47 Hz to 63 Hz
Number of treatments per charge	Approx. up to 8 treatments when nebulizing 2.5 ml of medication using a PARI Reusable Nebulizer	Power wattage	18 watts under normal load
Normal time to charge	Approx. 2 hours	Complies with	IEC 60601-1/IEC 60601-1-2
Ambient temperature range	50°F to 104°F (10°C to 40°C)		

**FEATURES OF THE PARI TREK® S COMPACT COMPRESSOR NEBULIZER SYSTEM**

# PARI TREK<sup>®</sup>S

COMPACT COMPRESSOR NEBULIZER SYSTEM

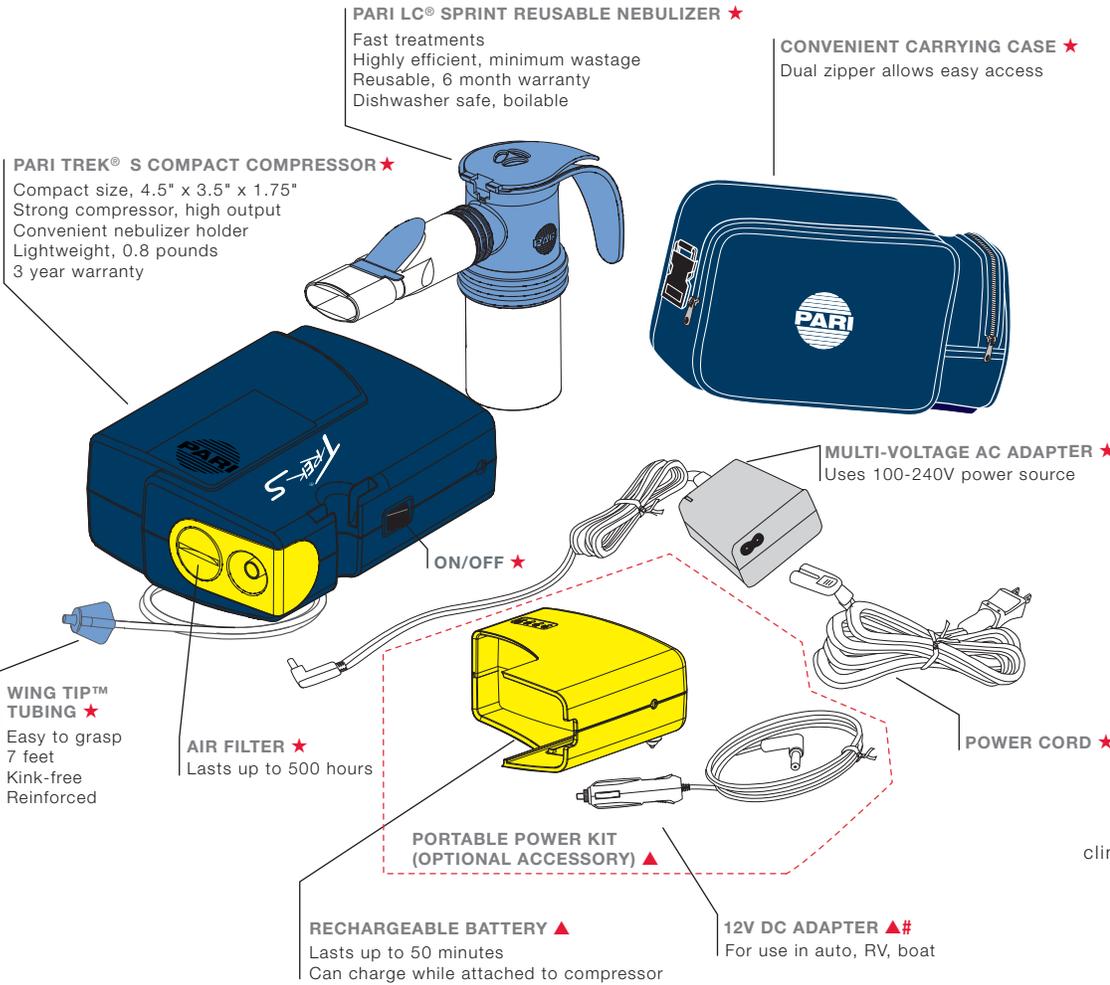

**PARI LC® SPRINT REUSABLE NEBULIZER ★**  
 Fast treatments  
 Highly efficient, minimum wastage  
 Reusable, 6 month warranty  
 Dishwasher safe, boilable

**CONVENIENT CARRYING CASE ★**  
 Dual zipper allows easy access

**PARI TREK® S COMPACT COMPRESSOR ★**  
 Compact size, 4.5" x 3.5" x 1.75"  
 Strong compressor, high output  
 Convenient nebulizer holder  
 Lightweight, 0.8 pounds  
 3 year warranty

**MULTI-VOLTAGE AC ADAPTER ★**  
 Uses 100-240V power source

**WING TIP™ TUBING ★**  
 Easy to grasp  
 7 feet  
 Kink-free  
 Reinforced

**AIR FILTER ★**  
 Lasts up to 500 hours

**PORTABLE POWER KIT (OPTIONAL ACCESSORY) ▲**

**RECHARGEABLE BATTERY ▲**  
 Lasts up to 50 minutes  
 Can charge while attached to compressor

**12V DC ADAPTER ▲#**  
 For use in auto, RV, boat

**PARI TREK® S PRODUCTS**

PARI TREK® S Combination Pack Part No. 47F35-LCS includes ▲ & ★ items

PARI TREK® S Compact Compressor with 12V DC Adapter Part No. 47F45-LCS includes ★ & #

PARI TREK® S Compact Compressor Part No. 47F10-LCS includes ★ items

PARI TREK® S Portable Power Kit Part No. 47F25 includes ▲ items

**PARI AND CLINICAL TRIALS**



PARI nebulizers are used in more clinical trials for nebulized medications than all other brands combined.


**PARI TREK® S PERFORMANCE WITH PARI LC® SPRINT REUSABLE NEBULIZER<sup>1</sup>**

<b>TOTAL OUTPUT RATE</b>	<b>.44 ml/min</b>
<b>MASS MEDIAN DIAMETER (MMD)</b>	<b>3.8 microns</b>
<b>% UNDER 5 MICRONS</b>	<b>65%</b>
<b>% UNDER 2 MICRONS</b>	<b>29%</b>
<b>AVERAGE TREATMENT TIME</b>	<b>5-6 min.</b>

<sup>1</sup>Measured with Malvern MasterSizer X at 50% relative humidity, 0.9% NaCl solution, inspiratory flow 20 liters/minute, continuous nebulization, 23°C, fill volume 2.5 ml.

**COMPARISON OF AEROSOLIZED MEDICINE DELIVERY DURING TYPICAL PATIENT BREATHING CYCLE**

Nebulized medicine is delivered more quickly into the patient's lungs with the breath enhanced PARI LC® SPRINT Reusable Nebulizer than with traditional nebulizers.

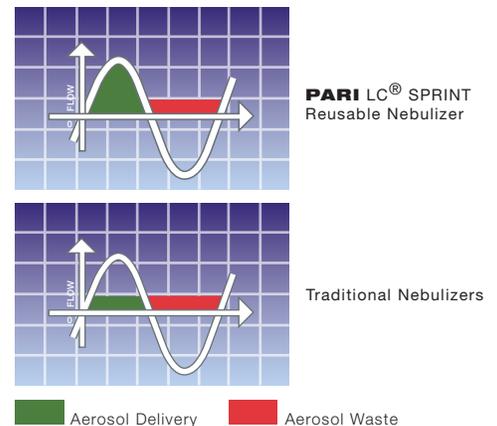

PART NO. 14F20-PVC MASK, ADULT (RIGID) PART NO. 14F20-PP (not shown)



BUBBLES THE FISH II™ PEDIATRIC AEROSOL MASK PART NO. 44F7248



PARI BABY™ CONVERSION KIT SIZE 1, 1-12 MTHS PART #44F2401 SIZE 2, 1-3 YRS PART #44F2402 SIZE 3, 3+ YRS PART #44F3301



Please register your product with PARI at [www.PARI.com](http://www.PARI.com)

[www.PARI.com](http://www.PARI.com) PARI RESPIRATORY EQUIPMENT, INC.  
 Email [ProductInfo@PARI.com](mailto:ProductInfo@PARI.com) 2943 Oak Lake Blvd. Midlothian, VA 23112-3998 USA  
 Respiratory CEU Programs Available USA and CANADA 1.800.FAST.NEB (327.8632)  
 © 2006 PARI Respiratory Equipment, Inc. Tel 1.804.253.PARI (7274) Fax 1.800.727.4112

